



60 fps High Speed Video intraoral Camera w/ Tiny Camera Head & Slim Design
Able to avoid blur photo when snapshot; Crystal Image quality w/ long vision range (5-300 mm focus)



5 Glass layers Lens,



8 pcs White LEDs



8.5 mm Height
12 mm Width


The Second Capture

Capture

Freeze / LED Power



Embedded ASCII Code Control for Buttons & Auto LED OFF/ ON (option)

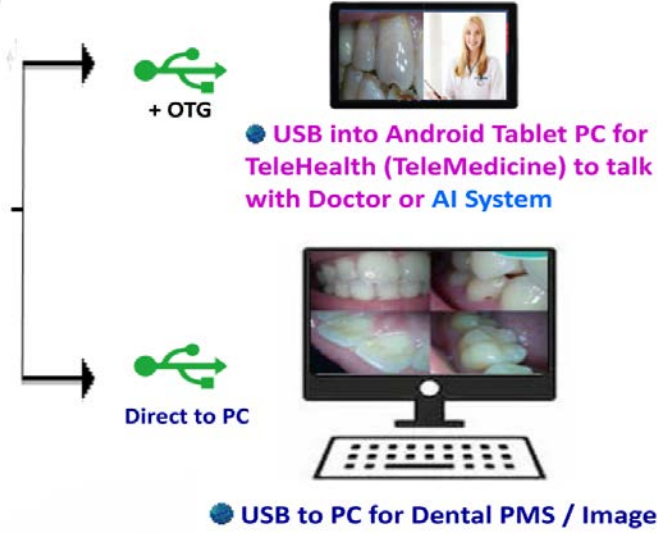

< Features >

- ◆ High Speed Video (60 fps) with more fluent video to reduce blur photo
- ◆ Tiny & Narrow Camera Head **12mm width & 8.5mm height**
- ◆ Five Glass Layers high-end Optical Camera Lens to work with Sony Sensor
- ◆ New Camera Lens to reach Super Image Quality w/ small Camera head
- ◆ Much less pressure when camera into patient's mouth
- ◆ Equipped with **ASCII code** to work with special Image Software
- ◆ A Convenient Camera Tool for **Telemedicine** or Telehealth
- ◆ Designed & Made in Taiwan
- ◆ 2.0 Mega Pixel sensor/ 1280 x 720 High Definition Video Streaming
- ◆ 110 ° View Angle & 5-300 mm long focal distance, a intra & Extra oral camera
- ◆ Good for one tooth, one row teeth, whole mouth and also Patient's face.
- ◆ "Freeze" function for Dentist to explain more to patients
- ◆ Auto LED Turn OFF / On when plugged IN / OUT from holder
- ◆ **USB 2.0** w/ long USB cable fixed into camera body

USB Cable fixed into Camera Body tightly, no more special Video connector; special connector is used to loose and high cost once changed